



## 2 FREQUENCIES = DOUBLE ACTION

To act simultaneously and selectively  
on 2 levels (superficial and deep fat).

ADVANCE **EB** BEAUTY



# ADVANCE BEAUTY

Rovaltain, parc 45 parallèle  
5 Rue Olivier de Serres  
26958 Valence Cedex 9  
FRANCE

Tél : +33 (0)4 75 48 17 54  
Fax : +33 (0)4 75 47 74 93

[Info@advance-beauty.com](mailto:Info@advance-beauty.com)  
[www.advance-beauty.com](http://www.advance-beauty.com)



Création : MACOMAVOI

MADE IN FRANCE  
BY ADVANCE BEAUTY

## THE NEW GENERATION **HAND FREE** SLIMING TREATMENTS



US2A **Sequential**



## TRANSDUCERS

The US2A Sequential comes with 12 ultrasound plates and each plate contains an assembly of 3 ceramic transducers. Therefore, with up to 36 points of impact on fat tissues, the US2A Sequential can treat the whole body. Low frequency ultrasound is used to reduce the ungraceful appearance of orange-peel skin.

The ceramic transducers used in US2A Sequential are specially made for our plates, thus they attain much better accuracy and efficiency than usual piezo electric cells ,as well as even power distribution to all areas of the plate.

This EXCLUSIVE advanced technology guarantee a perfect performance.  
Plates are associated in pack of 2 or 3.

## ULTRASOUND

Ultrasound works by mechanical and stimulative actions depending on the frequency and the emission power level. Choice of frequency delineates the depth of treatment. Low frequency ultrasound (USLF) at right intensity has a mechanical action on superficial fat. Medium frequency ultrasound (USMF) stimulates the sympathetic nervous system that induces lipolysis of the internal and subcutaneous fat layers.

The sonophoresis process, that happens when ultrasound is applied, activates the penetration of our Dermalyse Gel. Its active cosmetic ingredients act in synergy to increase slimming effectiveness.

### MECHANICAL & INFORMATIVE ACTIONS

- Activates the lipolysis process in only 30 mins
- Soften tissues and break collagen columns
- Large area treated at once
- Automatic treatment (Hand free device)
- Over 1 sqm treated at once
- Easy, comfortable and painless

### EXCLUSIVE CONCEPT

Sequential US2A comes with:


- 2 packs of 3 USLF plates
  - 2 packs of 2 USLF plates
- Action on 30 points of impact  
Superficial fat  
Clear the orange peel texture

- 1 pack of 2 USMF plates
- 6 emissions reactivating the neuro-sympathic system  
destocking of deep greases

### INNOVATION: SEQUENTIAL PATTERN

In case of stress (physical reaction coming from external stimuli) the natural processes managing the skin will be reactivated, improving micro-circulations and relaunching lipolysis process. In a classic use of ultrasound on the skin the stress process is very important during the first few minutes and through habituation loses its stimulating effects as the treatment progresses.

The Sequential mode principle, patented by Advance Beauty, is based on a chain of ultrasonic sequences, different and non repetitive. So the state of stress remains constant throughout the treatment to increase the results and guarantee maximum efficiency for the 100% of the session.

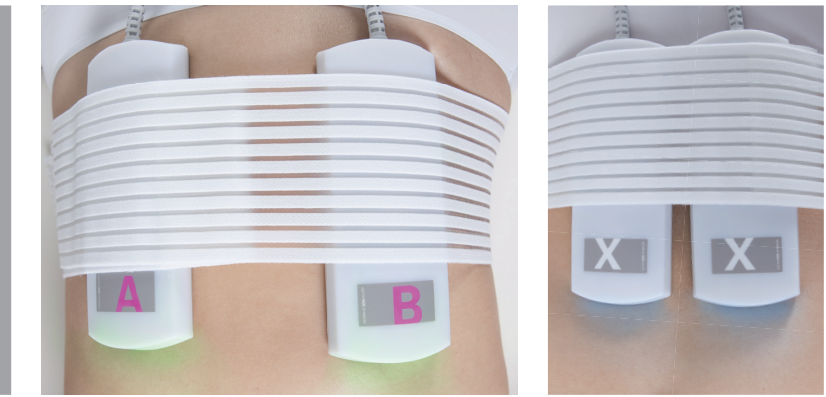

Design and user friendliness with interactive 9" screen.

- Only 5 cables are connected to the device
- Storage plates system for the packs

Outstanding profitability, thanks to its automatized «hand free» working for a session from 20 to 30 minutes



### TECHNICAL SPECIFICATIONS

Power supply	100-230 V AC / 50-60 Hz
Emission type	sequential, continuous or modulated
Maximum ultrasound power	2 W/cm <sup>2</sup>
Display 9" coloured	Screen with sensory touch
Programs	Intuitive sequences
Dimensions (cm)	25 (h) x 51 (L) x 30 (l)
Weight without plates	8 Kg
Easy storage plates system	Delivered (self assembly kit)
Connection to system	USB 2 type



### ACCESSORIES CHARACTERISTICS

Ultrasound L.F. plates	Nominal frequency	USLF - 30 KHz
	Number of cells / plate	3
	Number of pack of 3 plates	2
	Number of pack de 2 plates	2
	Total number of cells	30
	Active area / plate	110 cm <sup>2</sup> (17 x 6,5 cm)
Ultrasound M.F. plates	Operation indicator	integrated leds
	Nominal frequency	USMF - 500 KHz
	Number of cells / plate	3
	Number of pack de 2 plates	1
	Total number of cells	6
	Active area / plate	110 cm <sup>2</sup> (17 x 6,5 cm)
	Operation indicator	integrated leds

2 sets of elastic belts having various lengths are delivered with Sequential US2A to fix plates in place.

# DERMA-LYSE

## PARABEN FREE

The Derma-Lyse Gel is particularly recommended for daily application on areas with extra fat and presenting an orange-peel skin look. Only one product providing 3 major actions.

1

## LIPOLYSIS

The lipolysis corresponds to the hydrolysis of triglycerides into fatty acids and glycerol. The Derma-lyse Gel is a professional formula highly concentrated in Caffeine and theobromine molecules that have a well-known action on the activation of the deterioration of triglycerides.

2

## DRAINAGE

Fucus vesiculosus is an extract of seagrass that favors the elimination of water and wastes. For people concerns by fluid retention, Derma-lyse is more particularly recommended.

3

## CIRCULATION

Thanks to the Aescine well-known for its strong circulatory properties, the Derma-lyse Gel favours the cutaneous microcirculation.

